

Kowa



Technology for Life Science

KOWA APPLANATION
TONOMETER

HA-2


HS
Optikmaschinen



Convenient Item For Daily Clinical Practice



● Standard set of HA-2

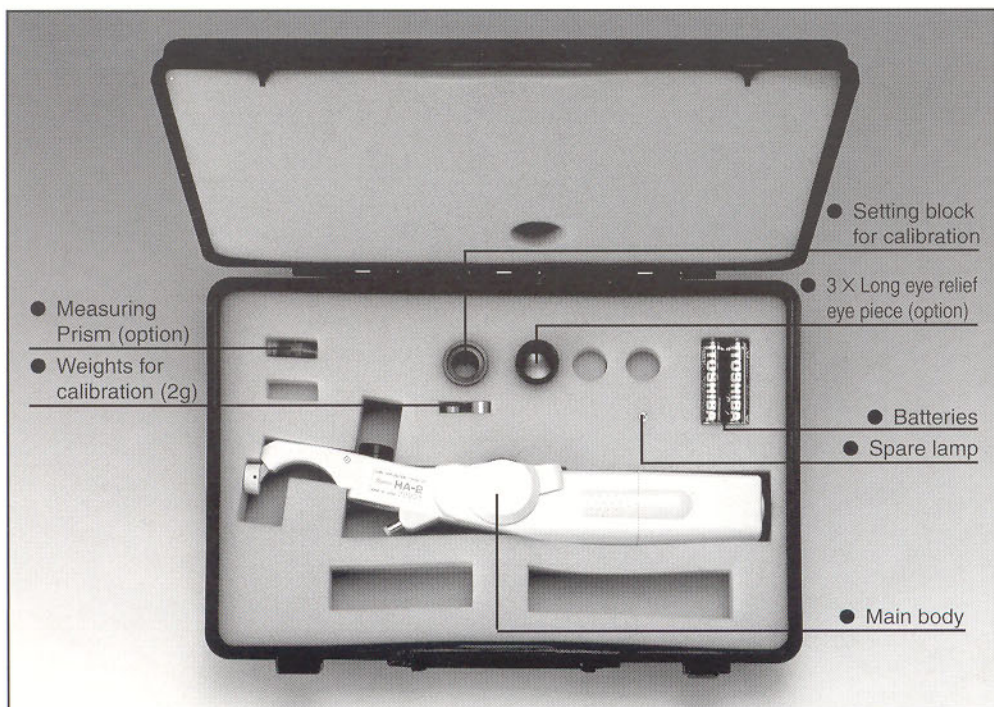

3 X Long eye relief eye piece (option)

● Main Features

- ◆ Compact and light-weight, well balanced and easy to handle by single hand.
- ◆ Unique one-spring mechanism assures correct reading regardless of patient posture.
- ◆ Built-in illumination permits you to use the instrument regardless of the illuminating condition of consultation room.
- ◆ Direct reading type precision dial requires no conversion table to change the readings.
- ◆ Accessories such as adjustable head rest and interchangeable eyepiece provide most suitable conditions for measurement.

● Specifications

Type	Hand-held applanation tonometer
Applanation prism	Doubling prism interchangeable (option)
Scale range	0~60 mm Hg (1 mm Hg division)
Illumination	
a) ON/OFF	Interlocked with thumb wheel
b) Bulb	Lens bulb 2.5 V, 0.25 A interchangeable
c) Battery used	1.5 V, 2 pcs Interchangeable
Eyepiece	
Magnification	6 X
Head rest	Movable
Length	290 mm
Weight	240g with batteries

* All specifications and external design are subject to modification of change without prior notice.



HS Optikmaschinen

Handwerkerstraße 14

48720 Rosendahl-Holtwick

Tel: 02566/4720

Fax: 02566/1620

Email: hsoptikmaschinen@hotmail.com

www.hs-optikmaschinen.de