



LED

CHART PROJECTOR

▶ PROGRAMMABLE TCP-1000

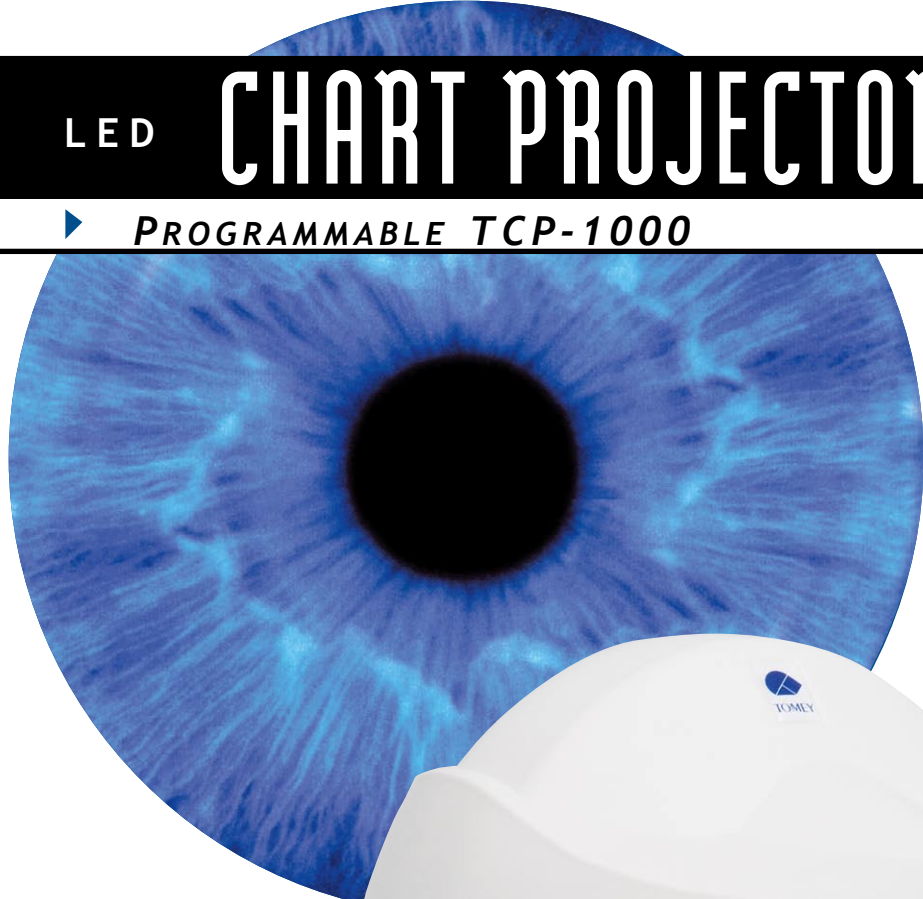


- *LED Light*
- *Brighter and Clearer Chart Display*
- *Sophisticated Design*
- *Quick Chart Selection*
- *Programmable Display Sequence*

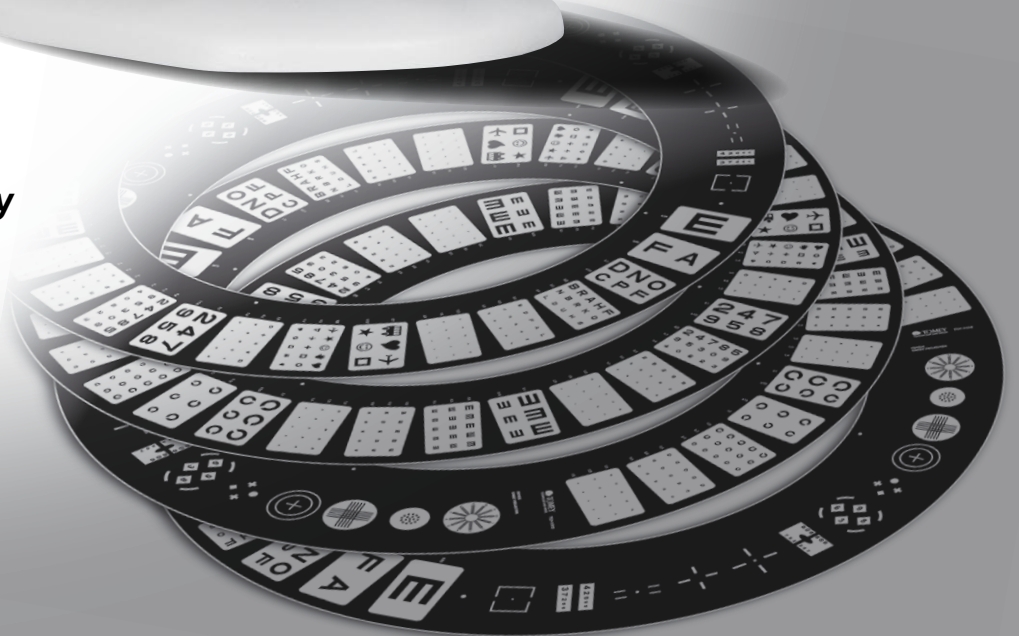

Modern Usability

The Tomey LED Chart Projector TCP-1000 is designed for durability and efficient operation.

Thanks to the LED light source the lamp is maintenance free due to its 500.000 hours lifetime. This makes the TCP-1000 very cost efficient and the chart displays are brighter and clearer because of the much better contrast. The intuitive remote control allows you to program three customised display sequences. Every sequence can be loaded with fifteen individual charts. You can choose them amongst 36 different test charts. The TCP-1000 will be delivered completely including remote control, projection screen and base for immediate use. You also can easily mount the system on your refraction unit or on the wall.

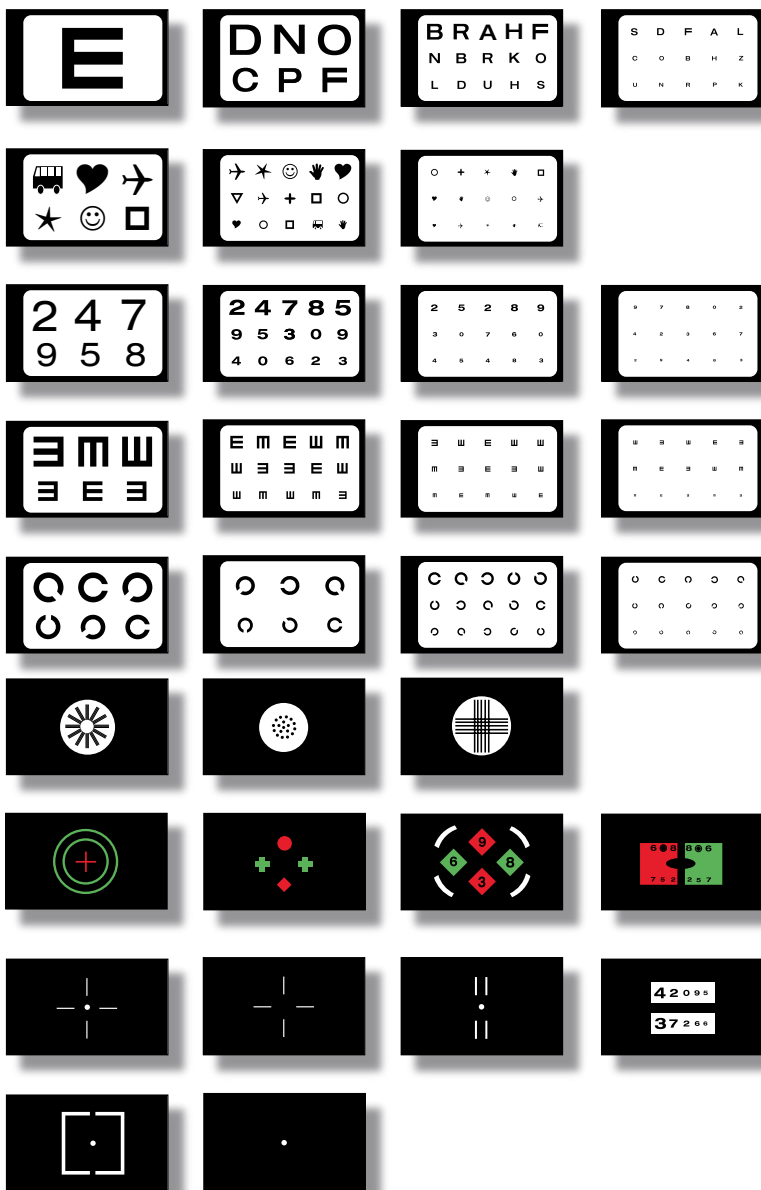

SPECIFICATIONS

Chart Change Time	0.014 s
Auto Off Time	5/10/15 min
3 User Programs	15 Steps Each
Masks	Single, Horizontal, Vertical, Red/Green, Polarizing
Projection Range	2.9 to 6.1 m
Projection Magnification	30x (at 5 m)
Lamp	3 W LED (373 to 801 nm)
Background Luminance	200 cd/m ² at 5 m (adjustable)
Projection Tilt	Approx ± 10°
Dimensions WDH	191 x 310 x 200 cm
Weight	4 kg
Power Consumption	11 W
Power Supply	AC 115V to 230V
Frequency	50/60Hz
Remote Control	Included (3 Channels Adjustable)